



time**to**relax

EYELOOP

BEDIENUNGSANLEITUNG | *MANUAL*



< relax in style

Since being established in 2000, the name CASADA stands for high quality products on the health-spa and fitness market. Today, CASADA-equipment is available in 37 countries world-wide.

Constant further development, an unmistakable design and a focused orientation towards a better quality of life are what characterize CASADA-products, as much as their unbeatable and reliable quality.

CONGRATULATIONS!

With the purchase of this massage device you have shown awareness of your health.

In order for you to enjoy the advantages of your massage device in the long term, we ask you to read and adhere to the instruction manual carefully.

We wish you a lot of fun with your personal EyeLoop.

17	Tradition and Modern
18	Safety instructions
19	Inside the package
19	Start up
19	Technical data
20	Equipment/functions
22	Remote control
24	Troubleshooting
24	Storage/care
24	Conformity declaration
25	Warranty conditions
26	Notes



The human eye is one of the most important sensory organs and it is a very reliable indicator of our overall physical well being. Nowadays physical comfort is often seriously impaired by strenuous computer tasks, professional stress, sleep and nutrient deficiency. One can quickly see these burdens in our eyes. They lose their expression and appear tired and exhausted.

The EyeLoop eye massager combines innovative technology and alternative Chinese acupressure. Based on this knowledge it gently relieves all blockages (red and sore eyes, swollen eyelids and headaches caused by tension) which come from stress and over exertion.

SAFETY INSTRUCTIONS



Please carefully read the following instructions to before using this fitness equipment to ensure its proper function and optimal effect. Please retain these instructions for use!

- This massage device complies with the recognised principles of technology and current safety regulations.
- This equipment is maintenance free. Any necessary repairs must be performed by an authorised professional. Improper use and unauthorised repairs are prohibited for safety reasons and will void the warranty.
- Never touch the power plug with wet hands.
- Please avoid the equipment coming into contact with water, high temperatures and direct sunlight.
- Be sure to unplug the equipment from the power point during extended periods of non-use to eliminate the risk of a short circuit. If you use the batteries: remove the batteries to prevent the leakage of battery fluid, which could cause errors.
- Do not use damaged plugs, cords or loose sockets. Damaged plug or cords must be replaced by the manufacturer, service representative or qualified personnel. In the event of malfunction unplug the equipment immediately. We assume no liability for damages resulting from wrongful or improper use.
- In order to prevent excessive stress on muscles and nerves, the recommended massage time of 30 minutes should not be exceeded.
- Never use pointy or sharp objects.
- Any massage – also massage by hand – has to be avoided during pregnancy and/or in case of one or more complaints in the area of the massage: Eye and head diseases recent injuries, thrombotic illnesses, inflammations and swelling of any kind as well as cancer. Consultation of your physician is recommended prior to any treatment of afflictions and ailments.
- Please remove contact lenses and hair jewellery before the massage.
- Only turn the device on, if you are already wearing it on your head to prevent damaged airbags by overpressure.
- Do not use the device while driving.
- If you depend on electrical aids such as a pace maker, please consult your **physician prior to usage** in any case.

INSIDE THE PACKAGE

- EyeLoop and remote control
- Adapter
- Batteries (3 x AA)

START UP

Connect the adapter to the EyeLoop and the power supply. For the mobile use: Remove the cover from the remote control and place the batteries in it. Close the cover again.

Tip: Please use either the adapter or battery. Do not use both at the same time.

TECHNICAL DATA

Measurements device: 22 x 10 x 8,6 cm

Measurements controller: 18,4 x 3,6 x 3,9 cm

Weight: 390 g

Rated power: 4,5 W

Battery: 3 x 1,5 V AA

Adapter: DC 4,5 V, 1000mA

Certificate:




EQUIPMENT / FUNCTIONS

We have numerous pressure points on the face and temples which in turn are connected to a respective area of the body, organ or muscle. With the help of an intelligent air-pressure control technology, Eye-Loop simulates a relaxing pressure massage with fingers, hands or palms exactly on these points. Its ergonomic design specifically adapts to your facial features and gives a relaxing eye and temple massage at all times. The gentle vibration stimulates the blood circulation, eases the feeling of tension and pressure on the eyes, reduces headaches and provides relaxing sleep.



AIR COMPRESSION

When using the Massage in the standard mode, the air pressure is reduced and the thermal function is set to a lower temperature. A massage in this mode helps the eyes to relax after a stresses and strains of everyday life. A three-stage interval mode brings an intensive massage experience. The mode change is done with a press of a button so that the pressure on the area being massaged is slowly increased. It alleviates the fatigue symptoms, relaxes and soothes the nerves and improves your general well-being.

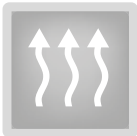


HEATING

The thermal function stimulates your metabolism and gives you a pleasant feeling during the massage.

The EyeLoop adapts perfectly to every head shape because of an adjustable strap. The simple regulation of the control panel allows for comfortable reprogramming without having to take of the massage goggles.

EQUIPMENT / FUNCTIONS

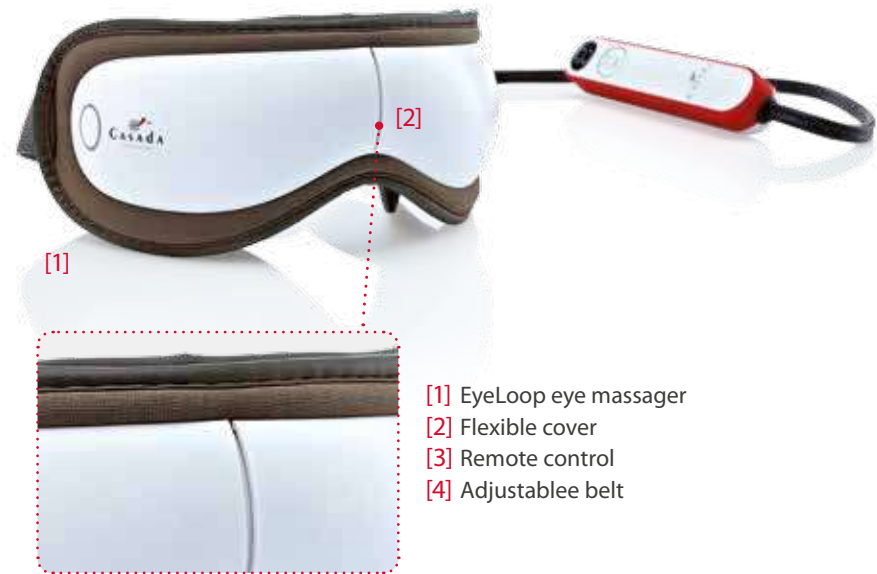

- [1] EyeLoop eye massager
- [2] Flexible cover
- [3] Remote control
- [4] Adjustable belt



- [4].....

- [3]

CASADA



- [1] Connection for power supply (for adapter)
- [2] Display for air pressure and heating mode
- [3] On/off
By using this button (holding down the button), EyeLoop starts. The default setting is the standard running time of 15 minutes, the air pressure mode 1 (relaxing) and the heating function in the lower temperatures (h). When the device is turned on this button can also adjust the warmth (see [4]). Hold the button down in order to turn off the EyeLoop.
- [4] On/off heating function
By using this button the heating function can be controlled among other things:
Press once: h = lower temperature
Press again: H = higher temperature
Pressing the button again will again display the selected air pressure mode.
- [5] Air Pressure
By using this button different air pressure modes can be selected. Select from three different modes with different intervals:
1: Relaxing
2: Kneading
3: Comfortable

TROUBLESHOOTING

Error	Adjustment
The massage device cannot operate even though it has been connected to the power supply.	<ul style="list-style-type: none">• Check the connection of the adapter.• Check whether the batteries are correctly installed.• The device could be overheated after prolonged use. It will automatically cool down to resume its normal condition. Please turn off the power and allow the device to cool off for 30 to 50 minutes.• Replace batteries.
Air compression is too weak.	<ul style="list-style-type: none">• Replace batteries.

STORAGE / CARE

Regular care will extend the life of your EyeLoop whilst also ensuring safe operation at all times.

1. Recommended storage: dry, free of dust.
2. Keep liquids and open fire away from the EyeLoop.
3. Clean the equipment using a dry, soft cloth.
4. Always unplug the EyeLoop when not in use.

CONFORMITY DECLARATION

Complies with European Regulatory Standards EMC 2014/30/EU and LVD 2014/35/EU.

WARRANTY CONDITIONS

The dealer grants the statutory warranty period on the product described overleaf. The warranty applies from the date of purchase. The date of purchase is demonstrated by the sales receipt.

The dealer is to repair or replace products registered within Germany which are found to be defective free of charge. Parts subject to wear and tear such as cloth covers are not covered. To assert your claim, simply return the defective product with the sales receipt before the end of the warranty period.

The warranty claim does not apply if it is established that a defect has occurred e.g. through external influence or as a result of repair or alteration not carried out by the manufacturer or an authorized dealer.

The warranty granted by the seller is limited to repair or replacement of the product. Under this warranty, the manufacturer or seller has no further liability and is not responsible for damage caused by disregard of the operating instructions and/or improper use of the product.

If inspection of the product by the seller reveals that the warranty claim relates to a defect not covered by the warranty or that the warranty period has expired, the cost of the inspection and repair are to be borne by the customer.



Casada International GmbH
Obermeiers Feld 3
33104 Paderborn, Germany
info@casada.com
www.casada.com

Copyright Images and texts are subject to the copyright of Casada International GmbH and may not be used for other purposes without express confirmation.

Copyright (c) 2017 Casada International GmbH. All rights reserved.

